

This Page Is Inserted by IFW Operations  
and is not a part of the Official Record

## **BEST AVAILABLE IMAGES**

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images may include (but are not limited to):

- BLACK BORDERS
- TEXT CUT OFF AT TOP, BOTTOM OR SIDES
- FADED TEXT
- ILLEGIBLE TEXT
- SKEWED/SLANTED IMAGES
- COLORED PHOTOS
- BLACK OR VERY BLACK AND WHITE DARK PHOTOS
- GRAY SCALE DOCUMENTS

**IMAGES ARE BEST AVAILABLE COPY.**

**As rescanning documents *will not* correct images,  
please do not report the images to the  
Image Problem Mailbox.**

**CADD-Diplomat PC Communications**

File Setup Help

Program Pump View Pump Reports View Saved Reports Outcome Data

Modem Control Disabled Direct Connect Disabled

Log-off

Check Program Check Accessories

### Communication Log

STOP START ENTER OPTIONS

LOCK NEXT Y

? DOSE N

**CADD Prizm**

Save Print

Not Connected Data Source: CADD Patient Data 3:43 PM

EXHIBIT

tabbles

B

**CADD-Diplomat PC Communications**

File Setup Help

Program Pump View Pump Reports View Saved Reports Outcome Data

Report Options:

Reports:

View

Print

Save

Not Connected Data Source: CADD Patient Data 3:42 PM

CADD-Diplomat PC Communications

File

Setup

Help

Program Pump

View Pump Reports

View Saved Reports

Outcome Data

View Options:

Open

Re-sort

Sorting Categories:

Reports will be sorted in the order listed below. To modify the sort order, select the desired category and use the UP and DOWN arrows to adjust the position.

↑

Move

↓

Patient

Report type

Report date

Pump serial number

Not Connected

Data Source: CADD Patient Data

3:44 PM

**CADD-Diplomat PC Communications**

File Setup Help

Program Pump View Pump Reports View Saved Reports **Outcome Data**

Name: ID:

No Patient Selected

Date of Birth: Phone Number:

Therapy: Start of Therapy: End of Therapy:

Complications: Diagnosis: Medication:

Access Device:

Patient Options:

Select


Add


Edit


Print

Delete


Not Connected Data Source: CADD Patient Data 3:45 PM

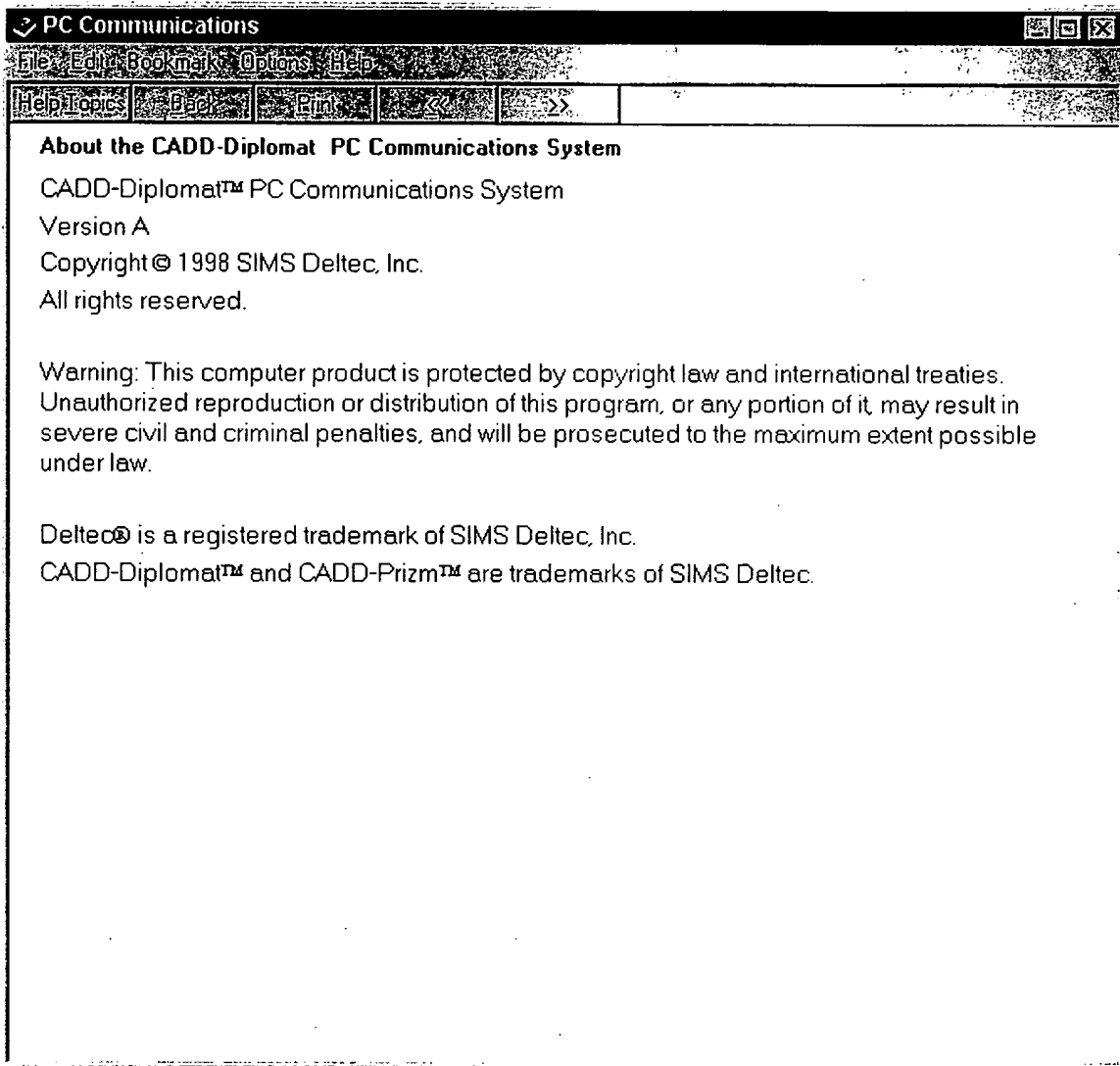


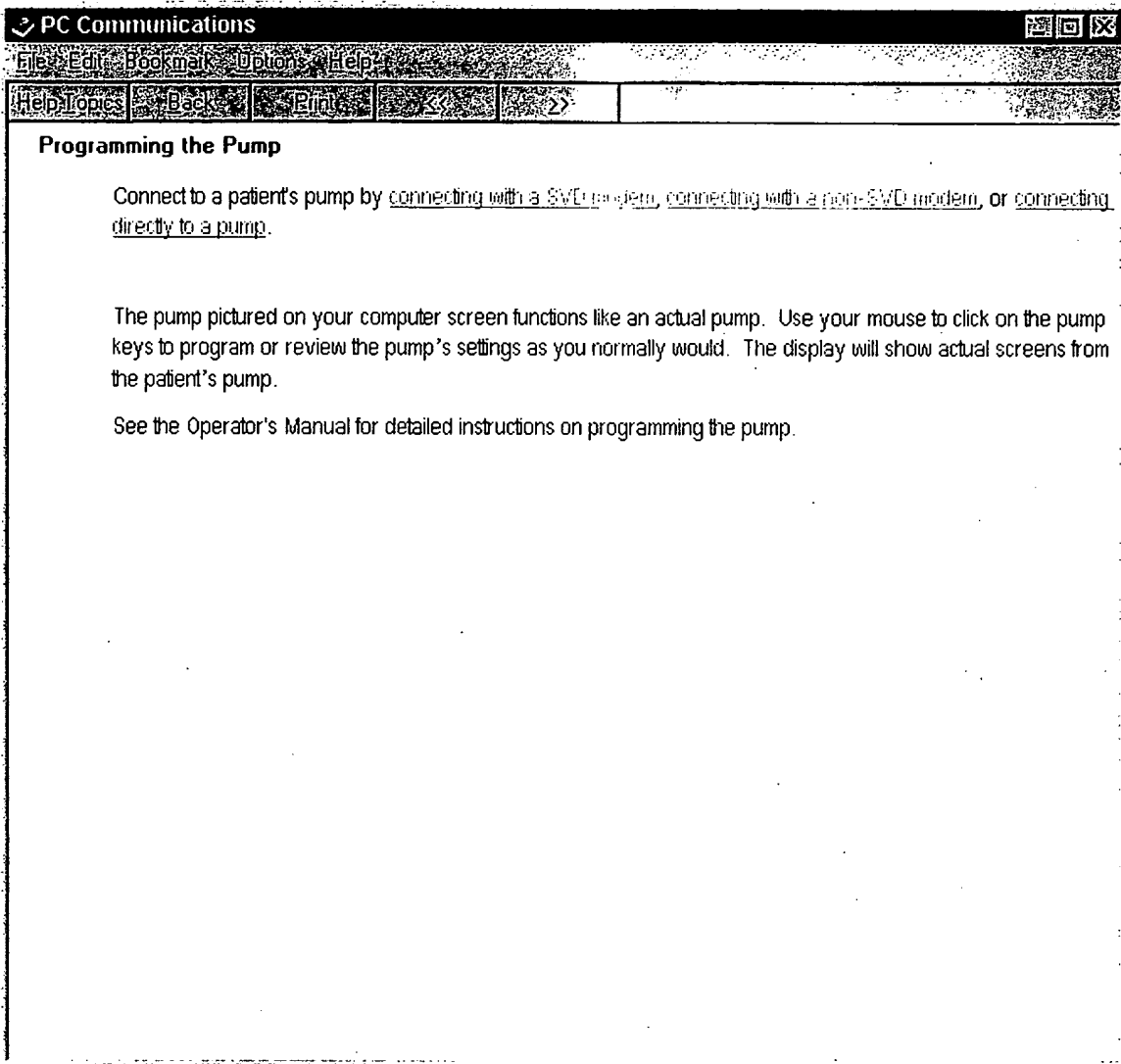

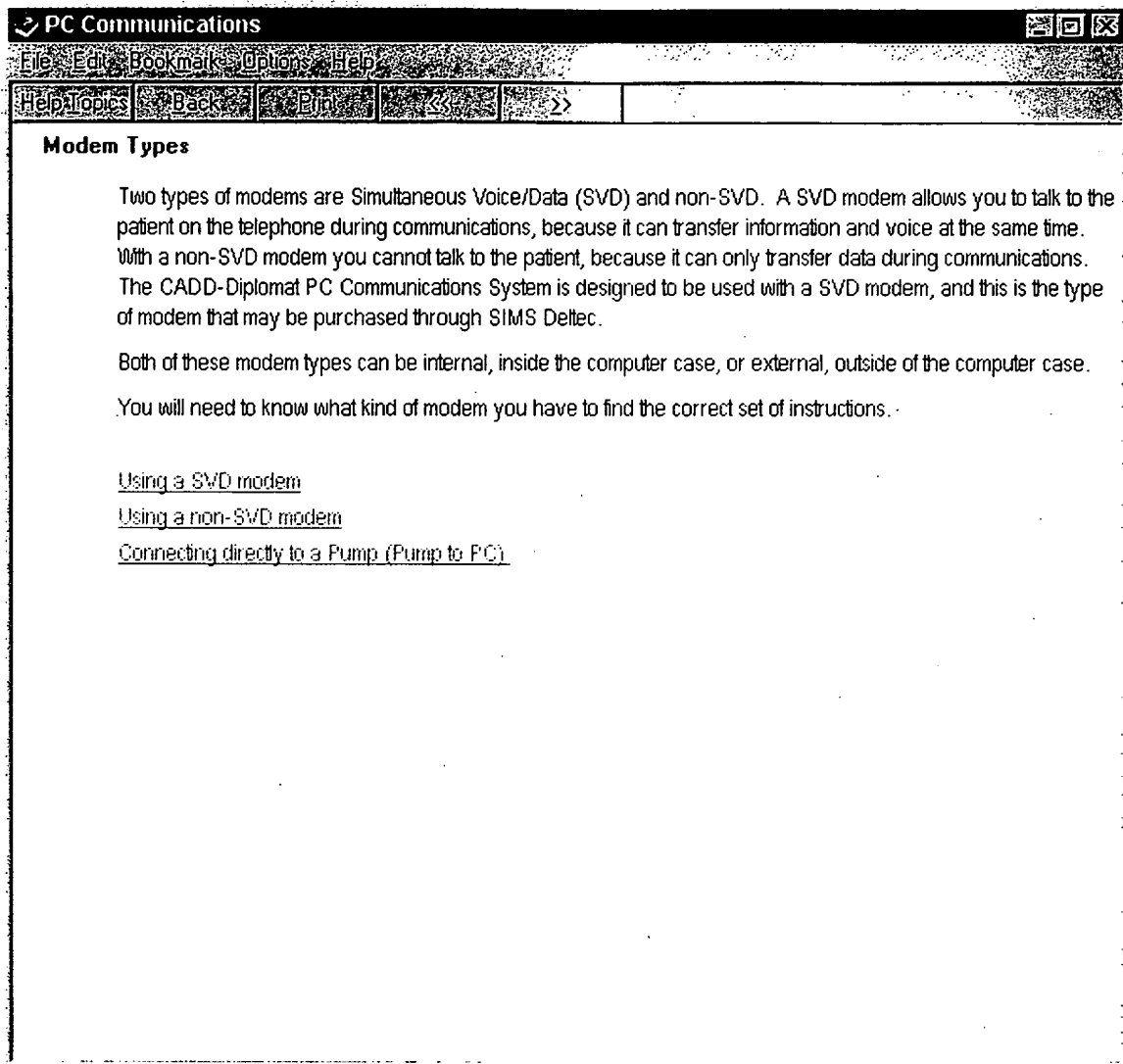



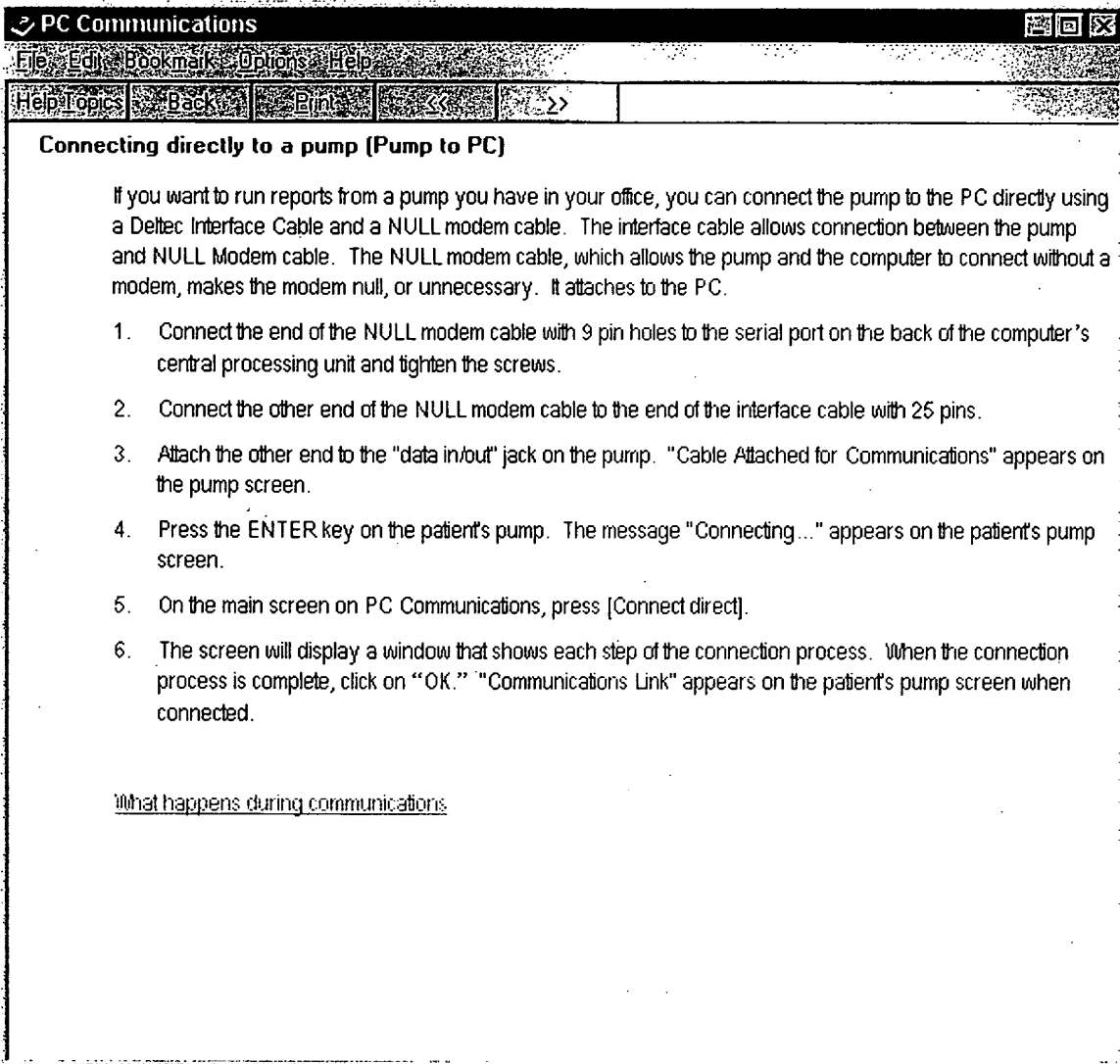

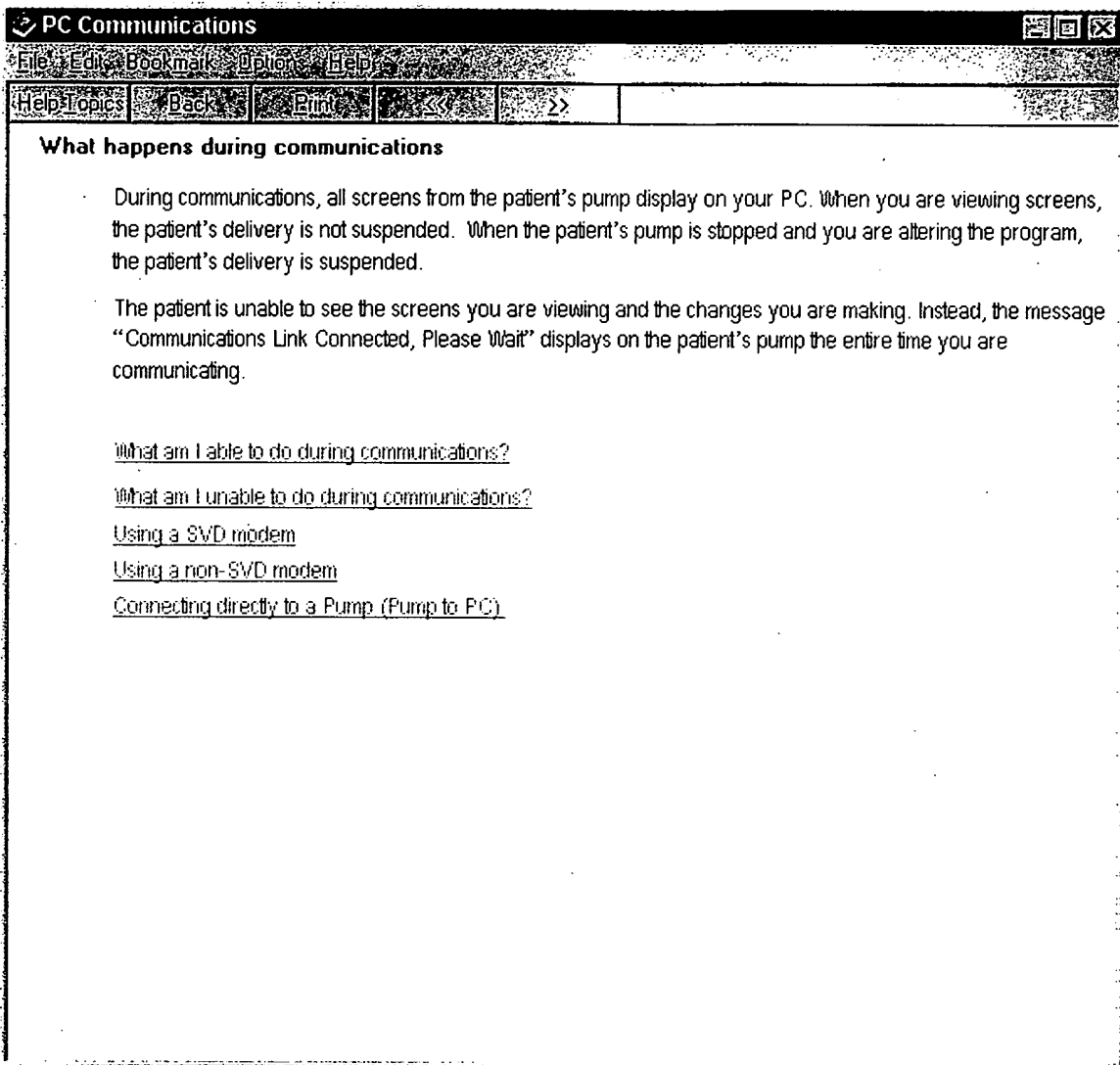

### What happens during communications

During communications, all screens from the patient's pump display on your PC. When you are viewing screens, the patient's delivery is not suspended. When the patient's pump is stopped and you are altering the program, the patient's delivery is suspended.

The patient is unable to see the screens you are viewing and the changes you are making. Instead, the message "Communications Link Connected, Please Wait" displays on the patient's pump the entire time you are communicating.


[What am I able to do during communications?](#)


[What am I unable to do during communications?](#)


[Using a SVD modem](#)

[Using a non-SVD modem](#)

[Connecting directly to a Pump \(Pump to PC\)](#)

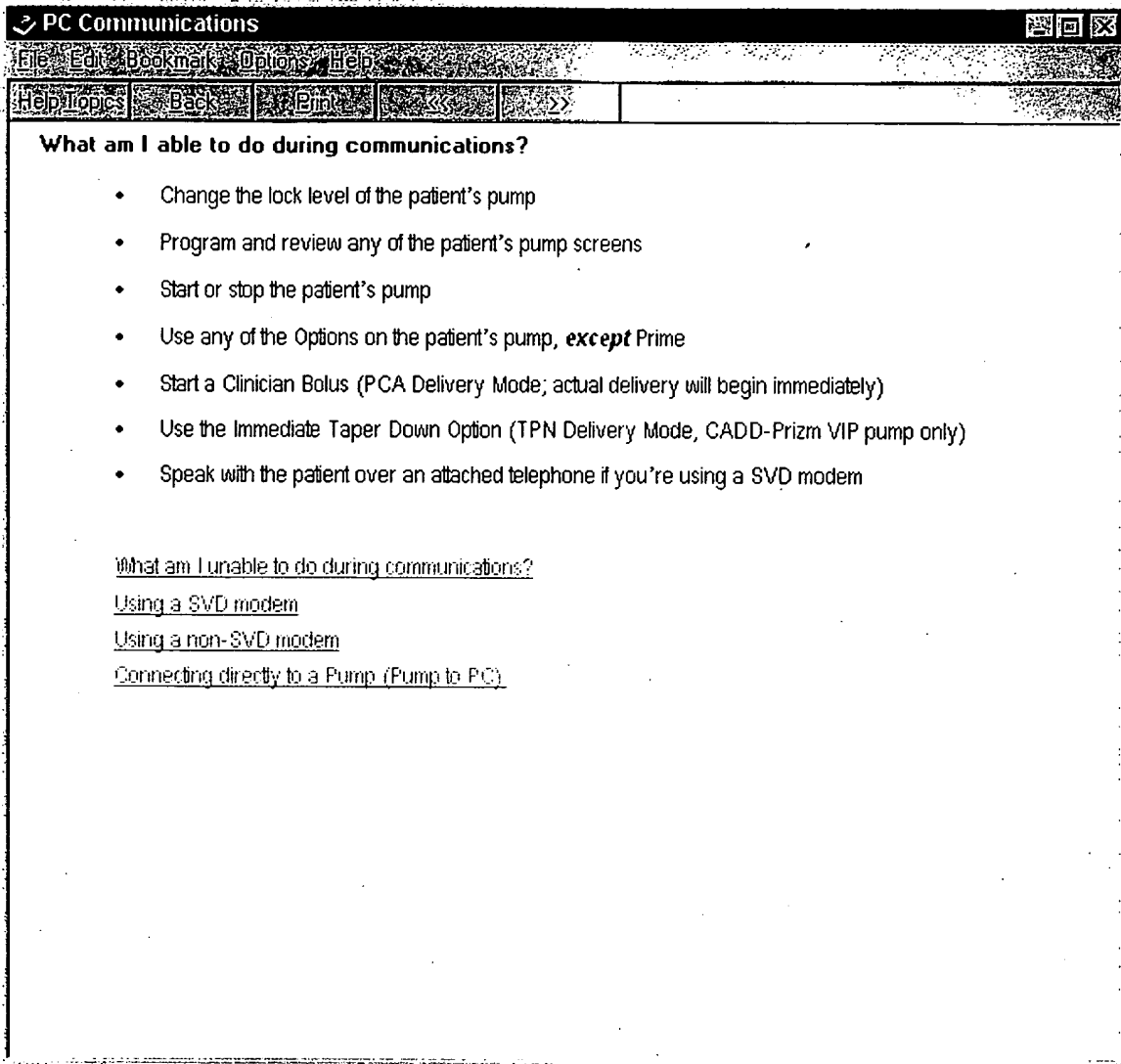

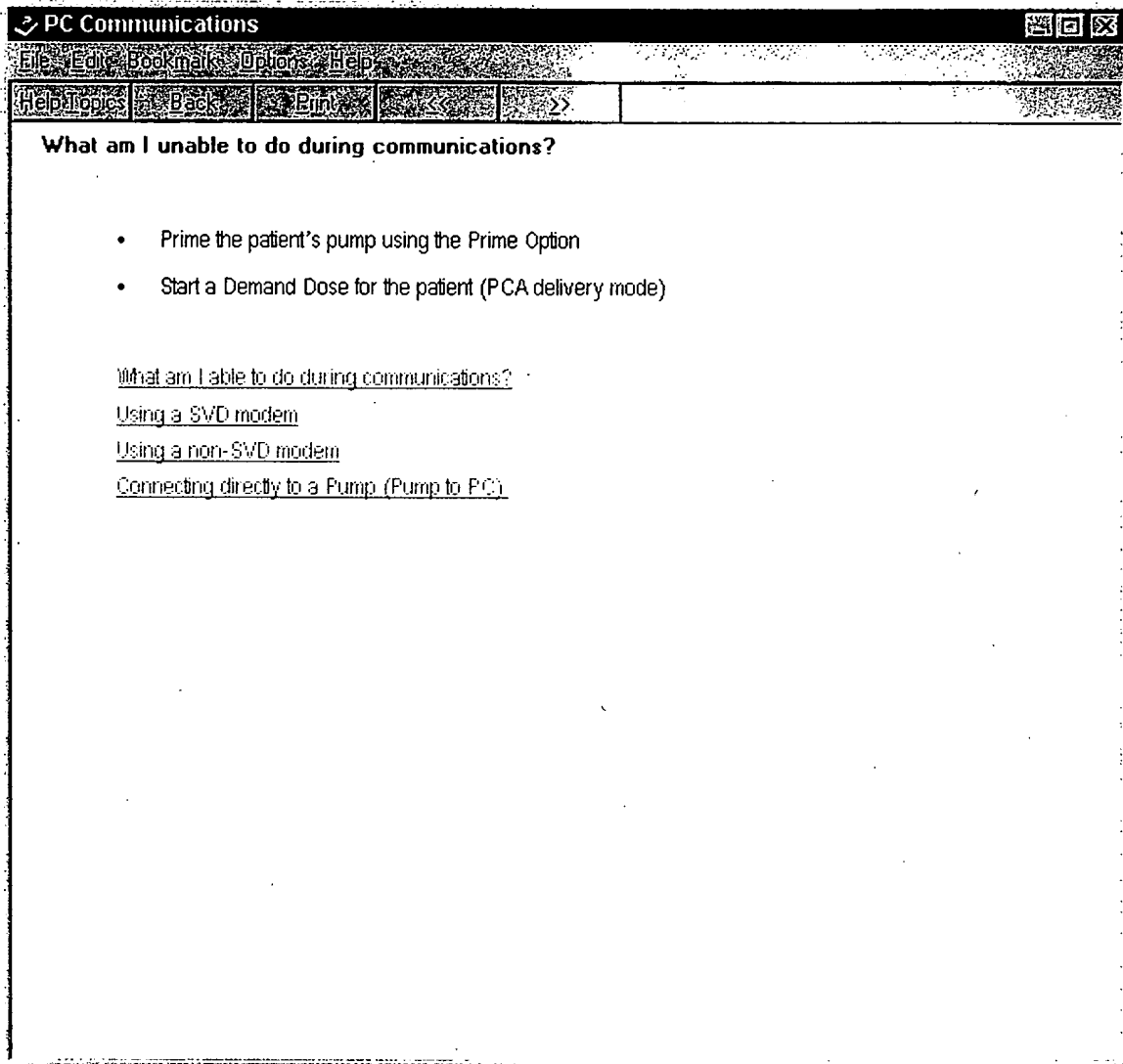

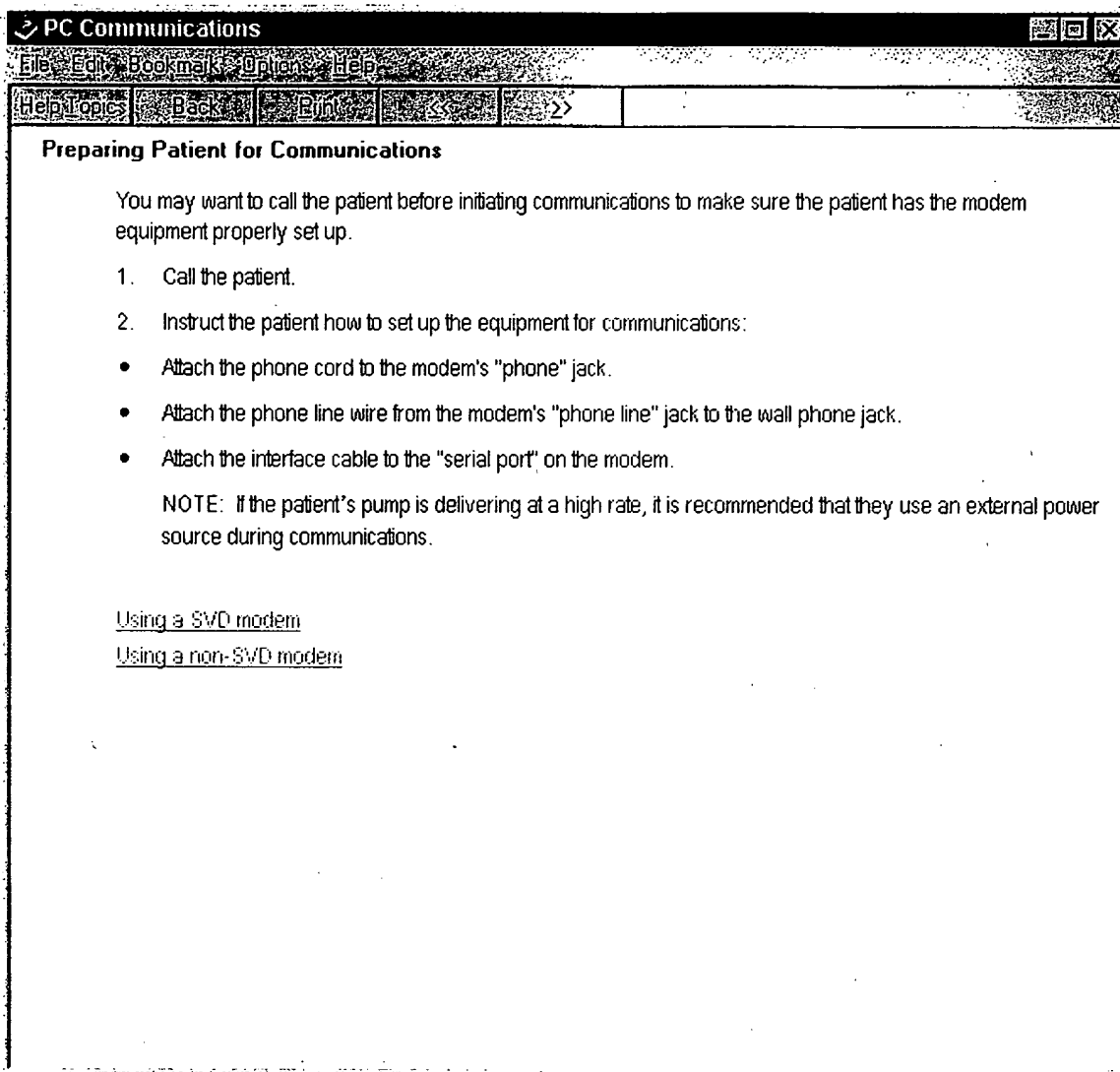












## Help Topics: PC Communications



Contents Index Find

Click a book, and then click Open. Or click another tab, such as Index.

-  **CADD-Diplomat PC Communications**
-  Password
-  System Setup
-  Disabling Call Waiting
-  Preparing Patient for Communications
-  Connecting to a patient's pump
-  Programming the Pump
-  Reports
-  Viewing or printing the communication log
-  Entering outcome data
-  Logging off from the pump

Open

Print...

Cancel

## Help Topics: PC Communications



Contents Index Find

1. Type the first few letters of the word you're looking for.

2. Click the index entry you want, and then click Display.

about the PC Communications System  
changing the password  
communications limitations  
communications overview  
connecting directly to a pump (Pump to PC)  
connecting with a SVD modem  
Diplomat limitations  
disabling call waiting  
entering outcome data  
important notes about PC Communications  
installing the modem  
logging off from the pump  
modem  
modem installation  
modem setup  
modem types  
non-SVD modem

Display

Print

Cancel

Help Topics: PC Communications



Content Index Find

1. Type the first few letters of the word you're looking for.

2. Click the index entry you want, and then click Display.

- modem types
- non-SVD modem
- notes about Diplomat
- overview
- overview of communications
- preparing patient for communications
- programming the pump
- removing the password
- serial port setup
- SVD modem connection
- using a non-SVD modem
- using help
- viewing or printing reports
- viewing or printing the communication log
- what am I able to do during communications
- what am I unable to do during communications
- working with saved reports

Display

Print...

Cancel