

This Page Is Inserted by IFW Operations  
and is not a part of the Official Record

## **BEST AVAILABLE IMAGES**

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images may include (but are not limited to):

- BLACK BORDERS
- TEXT CUT OFF AT TOP, BOTTOM OR SIDES
- FADED TEXT
- ILLEGIBLE TEXT
- SKEWED/SLANTED IMAGES
- COLORED PHOTOS
- BLACK OR VERY BLACK AND WHITE DARK PHOTOS
- GRAY SCALE DOCUMENTS

**IMAGES ARE BEST AVAILABLE COPY.**

**As rescanning documents *will not* correct images,  
please do not report the images to the  
Image Problem Mailbox.**

**CADD-Diplomat PC Communications**

File Setup Help

Program Pump View Pump Reports View Saved Reports Outcome Data

Modem Connect Disabled Direct Connect Disabled

Log-off

Check Program Check Accessories

### Communication Log

STOP  
ST/HL

ENTER

OPTIONS

LOCK


NEXT

Y

?

DOSE


N

CADD  Prizm

Save

Print

Not Connected Data Source: CADD Patient Data 3:43 PM

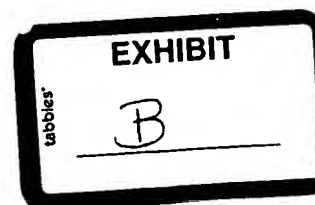

**CADD-Diplomat PC Communications**

File Setup Help

Program Pump View Pump Reports View Saved Reports Outcome Data

Report Options:

Reports:

View

Print

Save

Not Connected Data Source: CADD Patient Data 3:42 PM

CADD-Diplomat PC Communications

File

Setup

Help

Program Pump

View Pump Reports

View Saved Reports

Outcome Data

View Options:

Open

Re-sort

Sorting Categories:

Reports will be sorted in the order listed below. To modify the sort order, select the desired category and use the UP and DOWN arrows to adjust the position.

↑

Move

↓

Patient

Report type

Report date

Pump serial number

Not Connected

Data Source: CADD Patient Data

3:44 PM

**CADD-Diplomat PC Communications**

File Setup Help

Program Pump View Pump Reports View Saved Reports **Outcome Data**

Name: ID:

No Patient Selected

Date of Birth: Phone Number:

Therapy: Start of Therapy: End of Therapy:

Complications: Diagnosis: Medication:

Access Device:

Patient Options:

Select


Add


Edit


Print

Delete

Not Connected Data Source: CADD Patient Data 3:45 PM

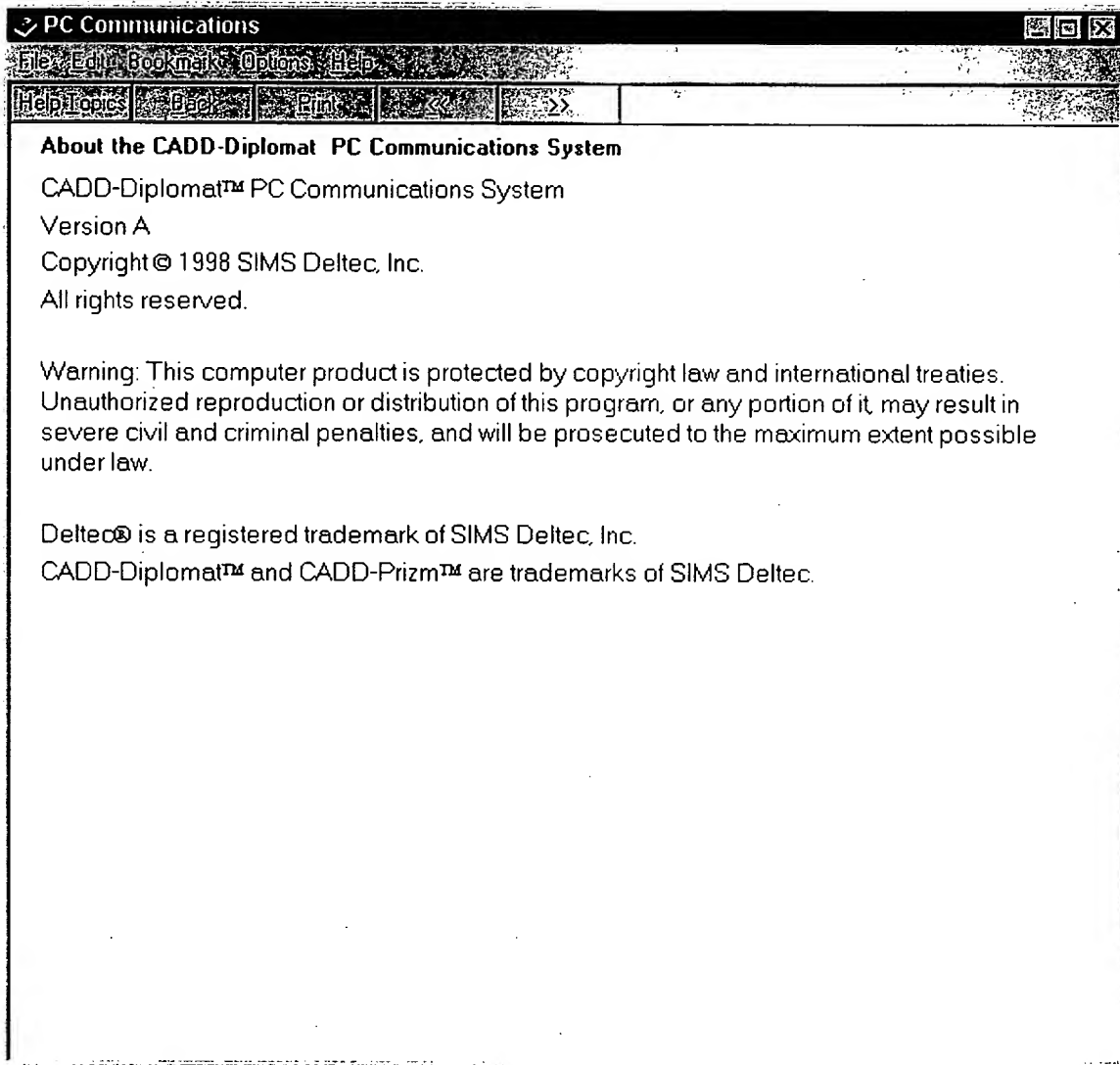

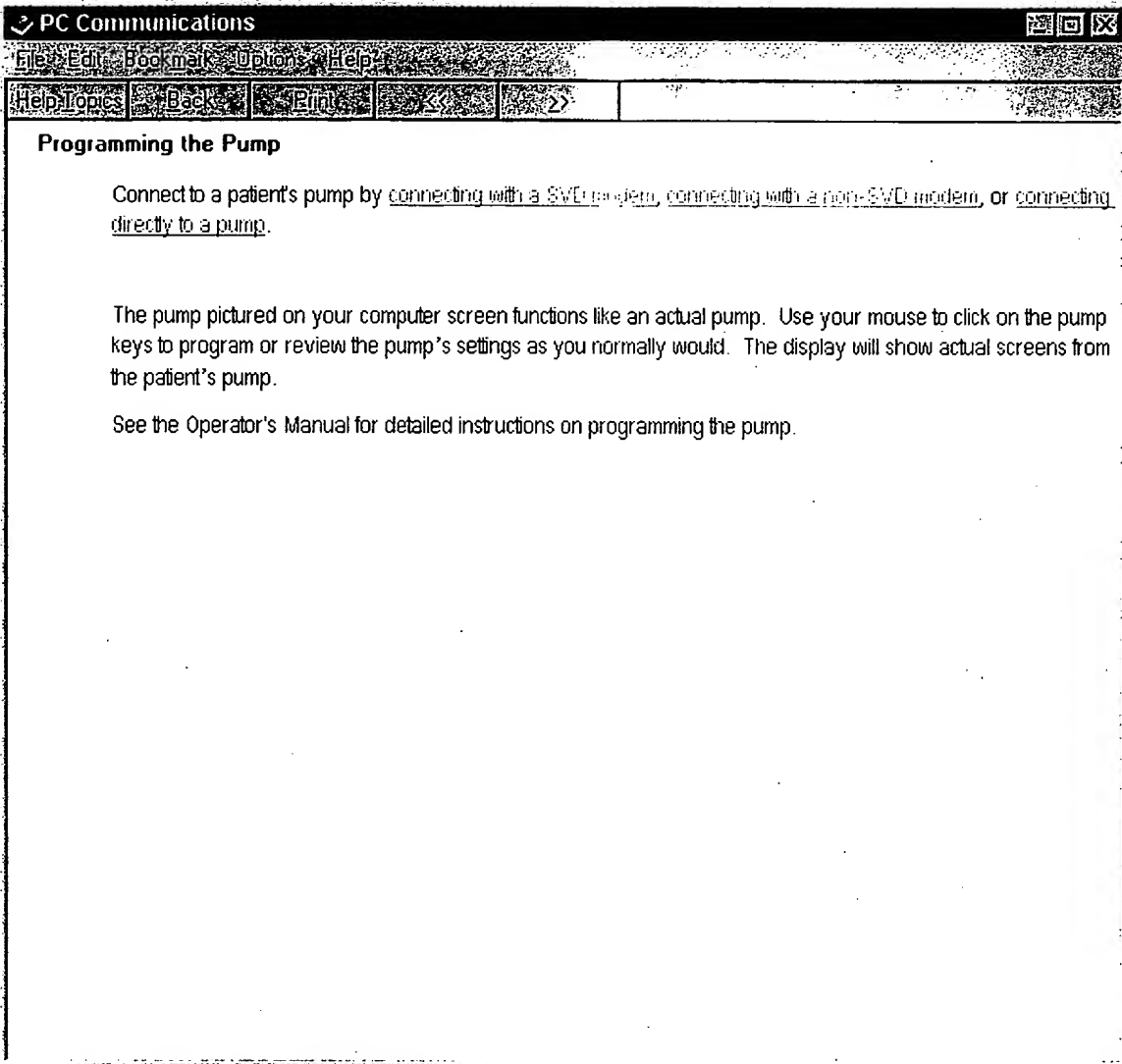

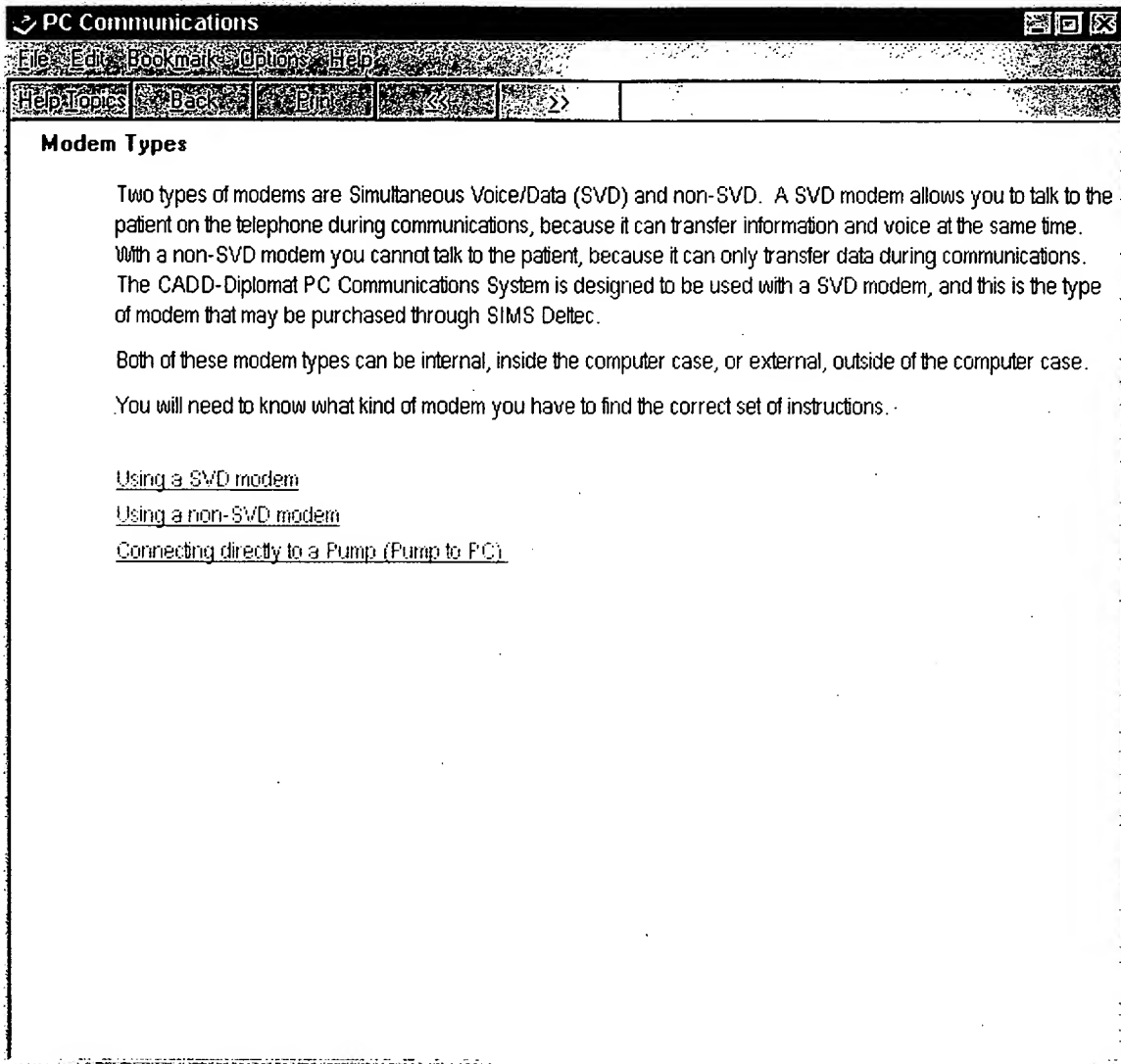

## **Modem Types**

Two types of modems are Simultaneous Voice/Data (SVD) and non-SVD. A SVD modem allows you to talk to the patient on the telephone during communications, because it can transfer information and voice at the same time. With a non-SVD modem you cannot talk to the patient, because it can only transfer data during communications. The CADD-Diplomat PC Communications System is designed to be used with a SVD modem, and this is the type of modem that may be purchased through SIMS Dettec.

Both of these modem types can be internal, inside the computer case, or external, outside of the computer case.


You will need to know what kind of modem you have to find the correct set of instructions.


[Using a SVD modem](#)

[Using a non-SVD modem](#)

[Connecting directly to a Pump \(Pump to PC\)](#)



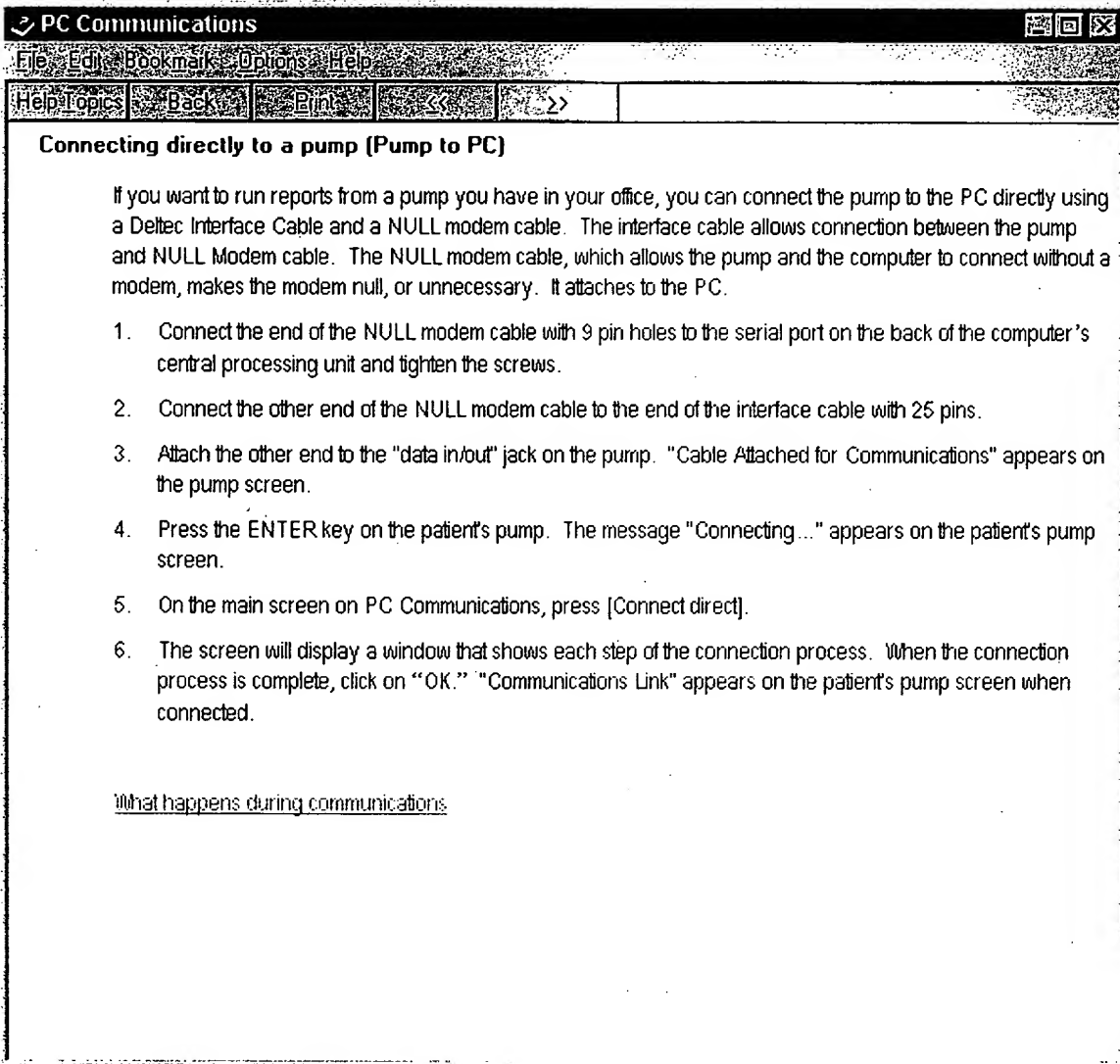

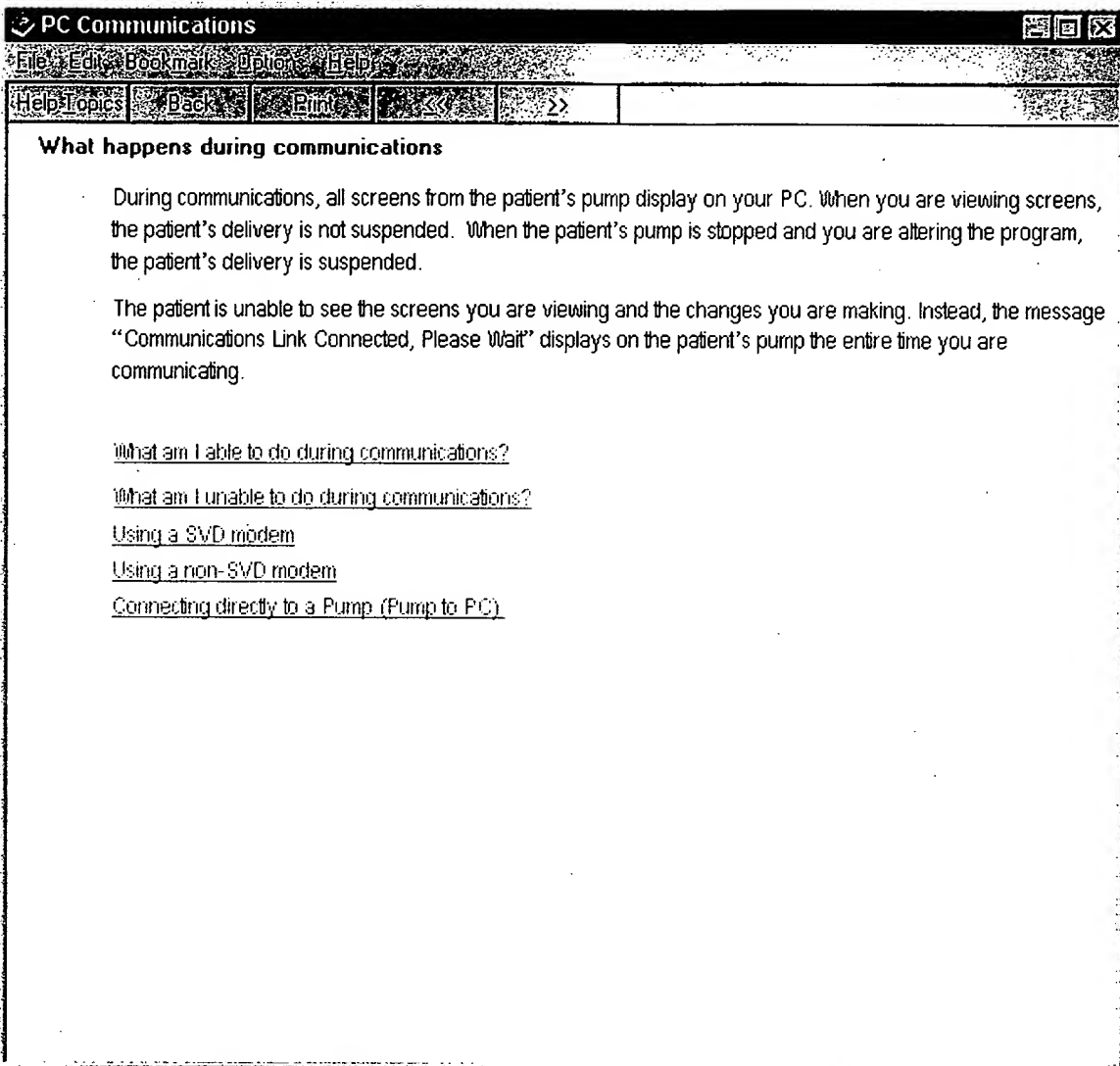

### What happens during communications

During communications, all screens from the patient's pump display on your PC. When you are viewing screens, the patient's delivery is not suspended. When the patient's pump is stopped and you are altering the program, the patient's delivery is suspended.

The patient is unable to see the screens you are viewing and the changes you are making. Instead, the message "Communications Link Connected, Please Wait" displays on the patient's pump the entire time you are communicating.


[What am I able to do during communications?](#)


[What am I unable to do during communications?](#)


[Using a SVD modem](#)

[Using a non-SVD modem](#)

[Connecting directly to a Pump \(Pump to PC\)](#)

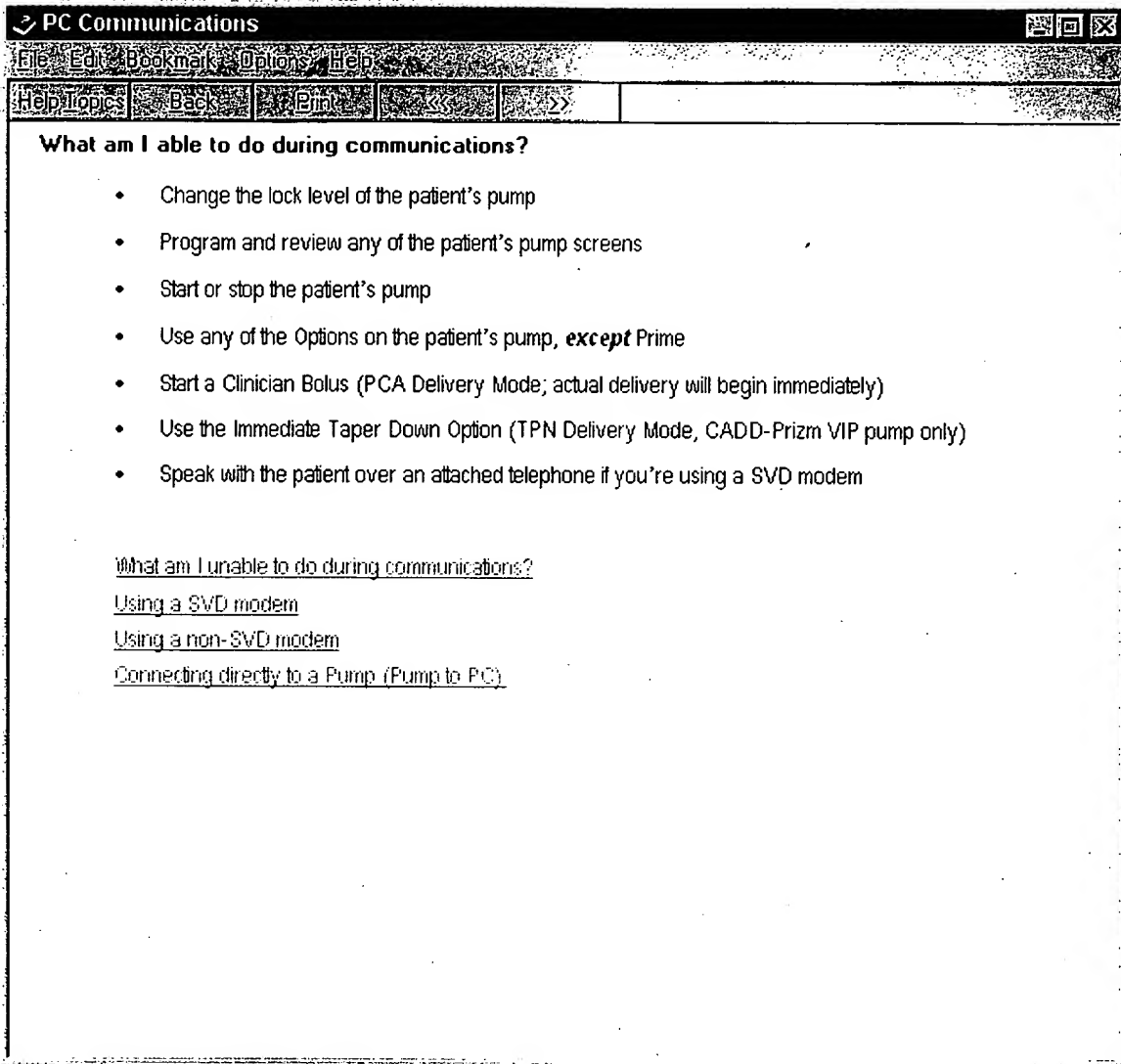

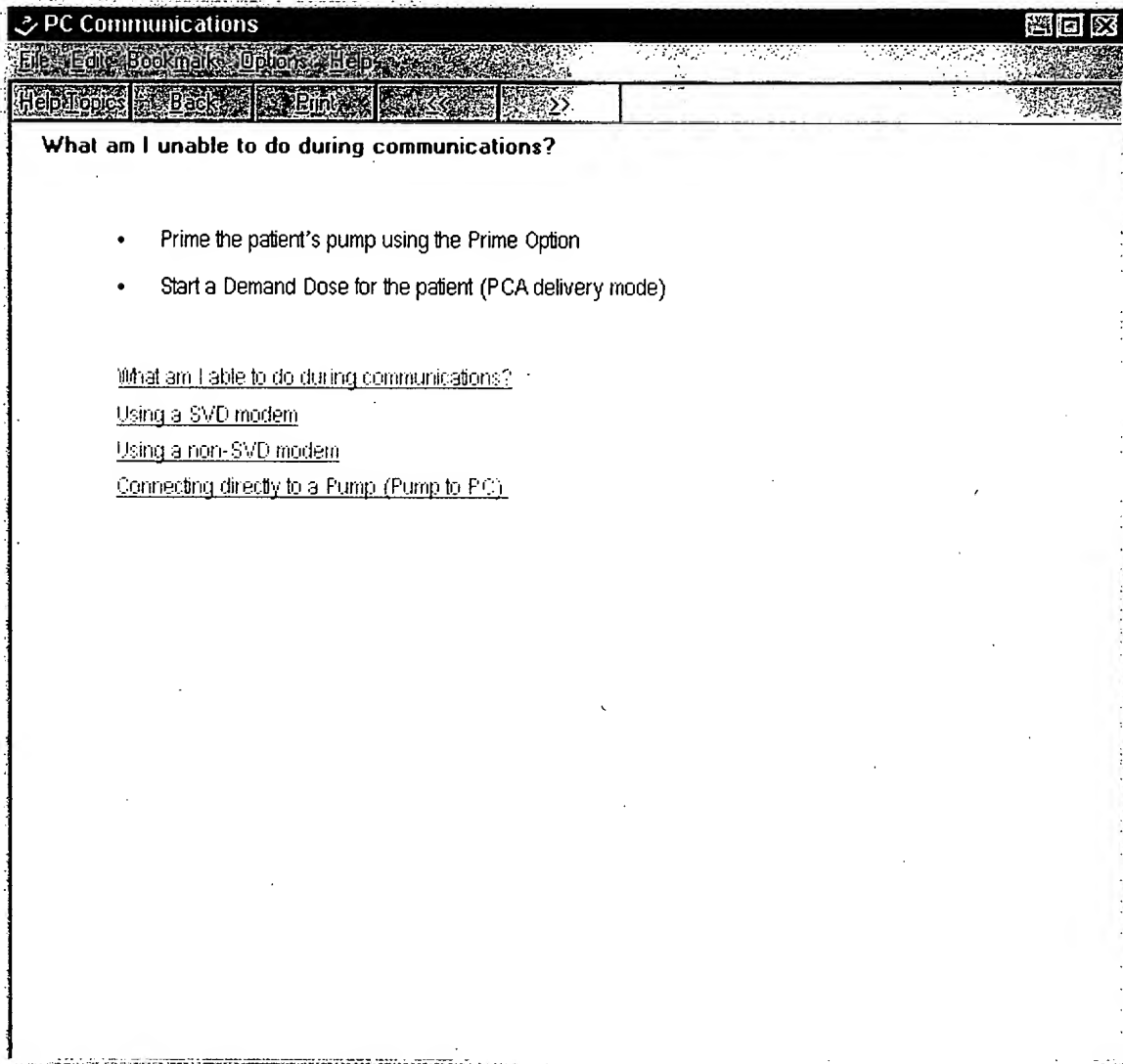

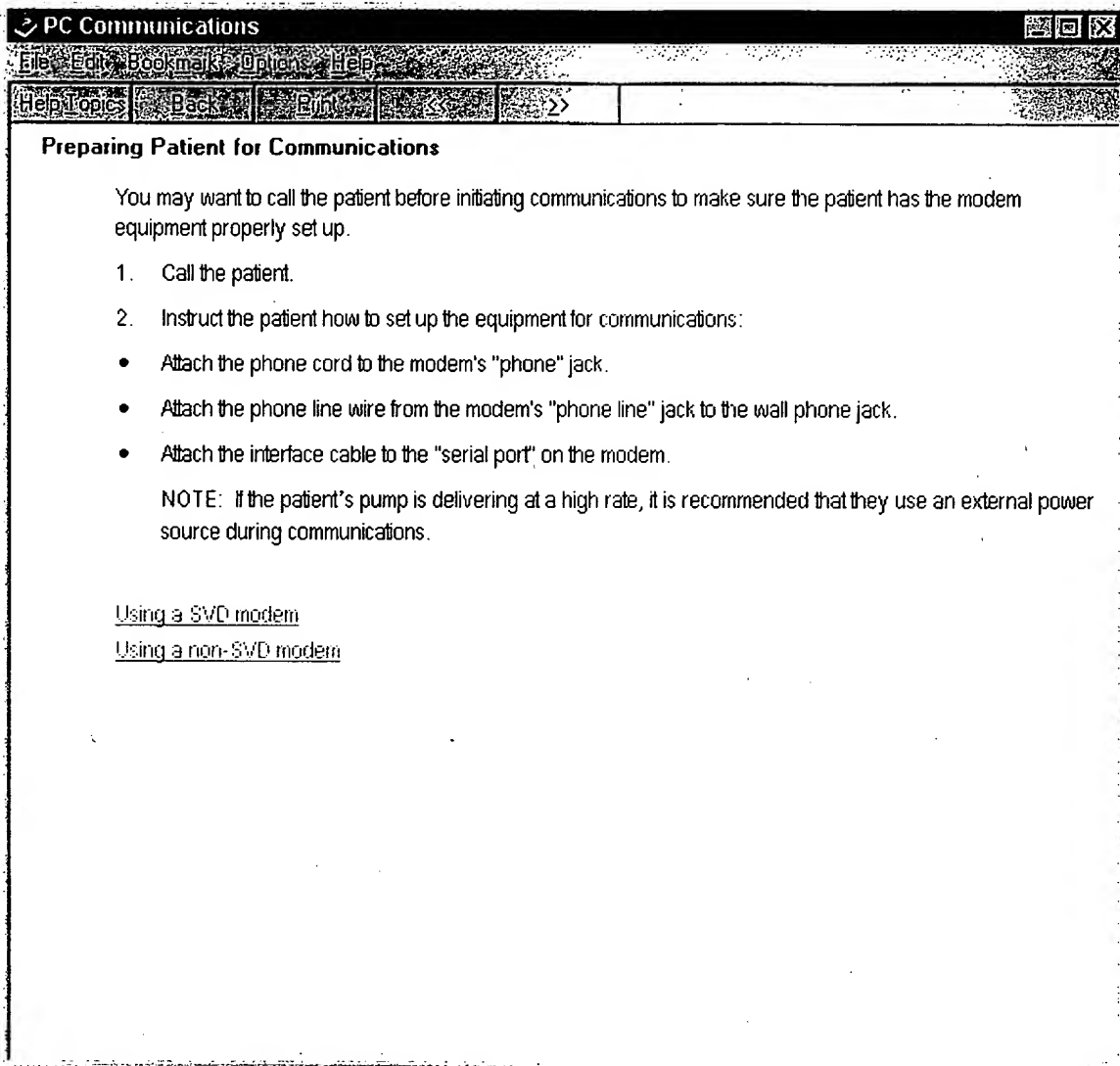

## Help Topics: PC Communications



Contents Index Find

Click a book, and then click Open. Or click another tab, such as Index.

- CADD-Diplomat PC Communications**
- Password
- System Setup
- Disabling Call Waiting
- Preparing Patient for Communications
- Connecting to a patient's pump
- Programming the Pump
- Reports
- Viewing or printing the communication log
- Entering outcome data
- Logging off from the pump

Open

Print...

Cancel

## Help Topics: PC Communications



Contents

Index

Find

1. Type the first few letters of the word you're looking for.

2. Click the index entry you want, and then click Display.

- about the PC Communications System
- changing the password
- communications limitations
- communications overview
- connecting directly to a pump (Pump to PC)
- connecting with a SVD modem
- Diplomat limitations
- disabling call waiting
- entering outcome data
- important notes about PC Communications
- installing the modem
- logging off from the pump
- modem
- modem installation
- modem setup
- modem types
- non-SVD modem

Display

Print

Cancel

Help Topics: PC Communications



Contents

Index

Find

1. Type the first few letters of the word you're looking for.

2. Click the index entry you want, and then click Display.

modem types  
non-SVD modem  
notes about Diplomat  
overview  
overview of communications  
preparing patient for communications  
programming the pump  
removing the password  
serial port setup  
SVD modem connection  
using a non-SVD modem  
using help  
viewing or printing reports  
viewing or printing the communication log  
what am I able to do during communications  
what am I unable to do during communications  
working with saved reports

Display

Print...

Cancel